

Ultrasonic

Homogenizer


20 KHz, 1200W



Homogenizer
Ultrasonic
1200W
20KHz

www.ultrasonic.co.ir

| <i>Property</i> | <i>UP1200</i> |
|---|-----------------|
| <i>Input Voltage (V)</i> | 220 |
| <i>Input Current (A)</i> | 6 |
| <i>Input Frequency (Hz)</i> | 50-60 |
| <i>Output power (W)</i> | 100-1200 |
| <i>Output Frequency (KHz)</i> | 20±1 |
| <i>Power Adjustable</i> | Yes |
| <i>Frequency Adjustable</i> | No |
| <i>Interface</i> | TFT LCD (touch) |
| <i>Working Mode</i> | Pulse |
| <i>Pulse ON (s)</i> | 1-20 |
| <i>Pulse OFF (s)</i> | 1-10 |
| <i>Total duration for a cycle (min)</i> | 30 |
| <i>Sound enclosure box</i> | Yes |
| <i>Thermometer</i> | Yes (PT-100) |
| <i>Temperature range (°C)</i> | 0-100 |
| <i>Horn Material</i> | Ti-6Al-4V |
| <i>Horn final diameter (mm)</i> | 30 |
| <i>Working volume (CC)</i> | 2000-5000 |

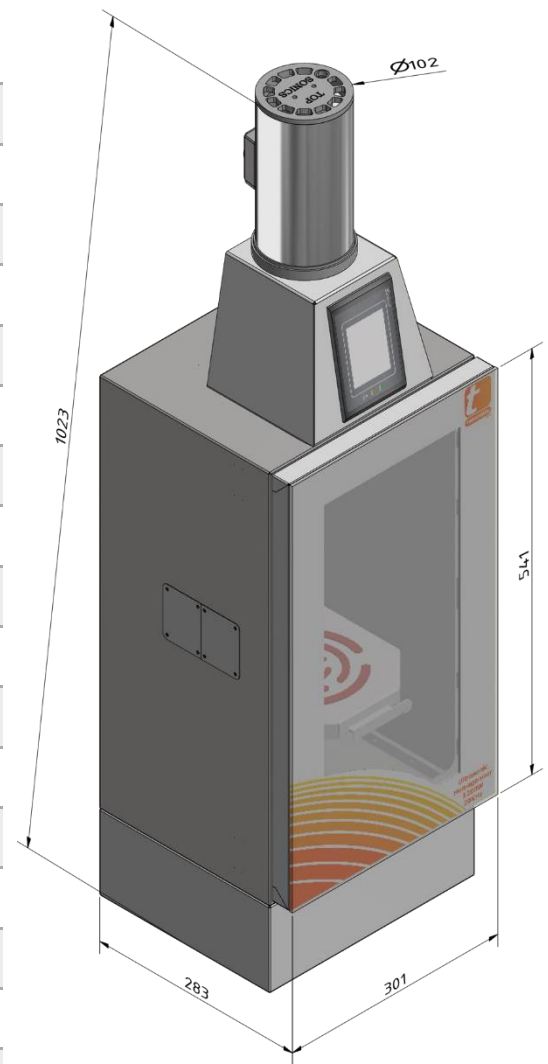

The scientific background of TOPSONICS dates back to 2003 with the design and manufacturing of ultrasonic surgical knife. Today, this company is one of the leading and fast growing in the field of power ultrasonic technology.

Ultrasonic waves are mechanical waves whose frequency exceeds human hearing range (20 Hz- 20 kHz.). Due to their property, these waves have different and somewhat interesting applications. Using a simple calculation, if a point with a frequency of 25 KHz. and a domain of 10 micrometer oscillates, its acceleration reaches a magnitude which is 25000 times greater than that of earth gravity ($g=9.8 \text{ m/s}^2$). Due to this acceleration and high speed, cavitation in liquids occur and during the phenomena of explosion of bubbles, an approximate pressure of 200 psi along with a temperature of 5000°K are observed.

On the other hand, if a relative movement with the above mentioned condition is occurred between two layers of thermoplastic, the high temperature leads to bonding of the two layers which is called ultrasonic welding one of the industrial application of ultrasonic technology.

Ultrasonic homogenizers are widely used in Nanotechnology to disperse nano material in to the oil and water. Besides that, phytochemistry, medicine, oil and gas are some other filed which ultrasonic homogenizer are useful.

آدرس: ایران، تهران، میدان توحید،
بیمارستان امام خمینی(ره)، مرکز رشد، اتاق
34، شرکت توسعه فناوری مافوق صوت.

00982166065887

00989123830551